Serie ACTIVE SCALES BLUE / CRYSTAL

MANUAL USUARIO / USER MANUAL / GUIDE DE L'UTILISATEUR /
MANUAL UTILIZADOR / MANUALE UTENTE / BEDIENUNGSANLEITUNG

ES / EN / FR / PT / IT / DE

BÁSCULA INTELIGENTE PARA LA GRASA CORPORAL

SMART BODY FAT SCALE

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EN

1. DOWNLOAD AND SET UP THE APPLICATION	10
2. BATTERY INSTALLATION	11
3. WEIGHT UNIT CONVERSION	11
4. CONNECT THE SCALE WITH YOUR PHONE	12
5. START MEASUREMENT	12
6. MULTIPLE USERS SUPPORTED	13
7. WARNING INDICATIONS	13
8. ABOUT THE BODY FAT SCALE	14
9. ADVISE FOR USE & CARE	14

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IT

1. DOWNLOAD AND SET UP THE APPLICATION	28
2. BATTERY INSTALLATION	29
3. WEIGHT UNIT CONVERSION	29
4. CONNECT THE SCALE WITH YOUR PHONE	30
5. START MEASUREMENT	30
6. MULTIPLE USERS SUPPORTED	31
7. WARNING INDICATIONS	31
8. ABOUT THE BODY FAT SCALE	32
9. ADVISE FOR USE & CARE	32

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1. DOWNLOAD AND SET UP THE APPLICATION

A. Go to the "App Store" or "Google Play", and search "SenssunLife" to download the Application.


B. Or in another most convenient way by scanning the QR code below.

SenssunLife


2. BATTERY INSTALLATION

Open the battery cover on the back of the scale as shown in the following picture, and then install the batteries according to the polarity mark in the battery case (Batteries non included).

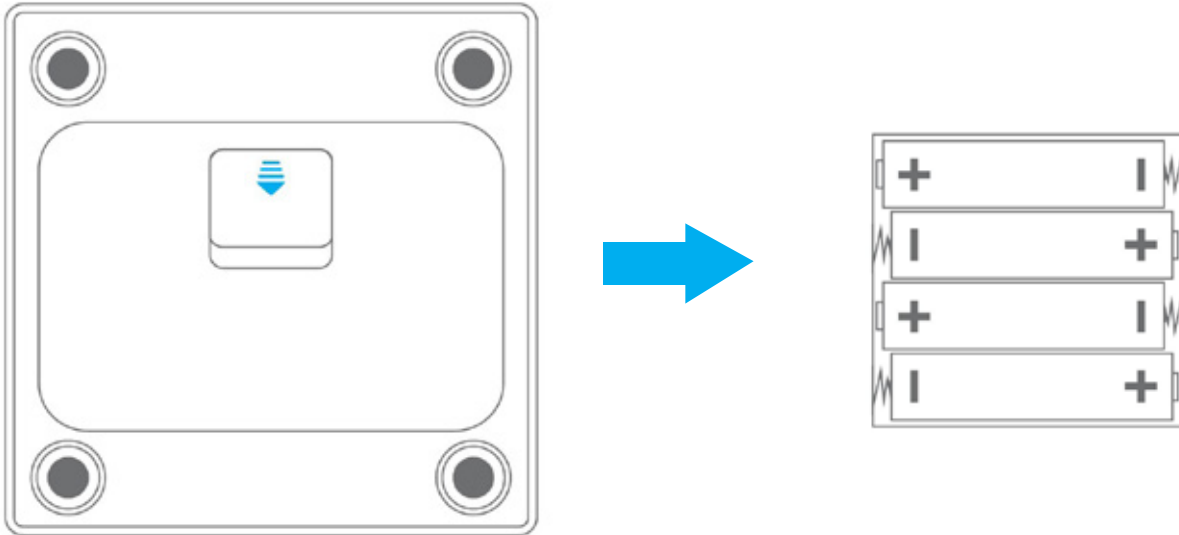3. WEIGHT UNIT CONVERSION


You will find the unit conversion button at the bottom of scale, press or pull it to convert the weight unit (Kg/lb), you will see it changed on LCD.

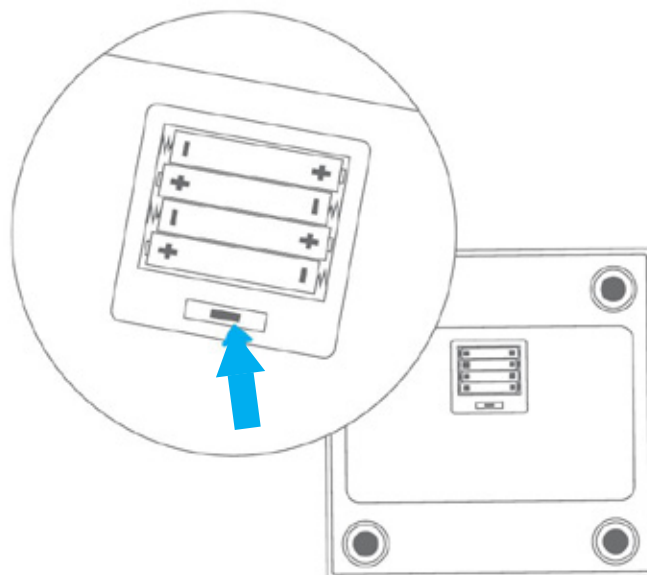4. CONNECT THE SCALE WITH YOUR PHONE

Turn on the Bluetooth of your phone and open the App. Sign up your account and follow the guideline to set up your profile and connection step by step. Being aware of it, the scale has to be on when you start to connect with it.

All measurements you taken can be saved in history for detail review anytime anywhere, by getting access to your own account.

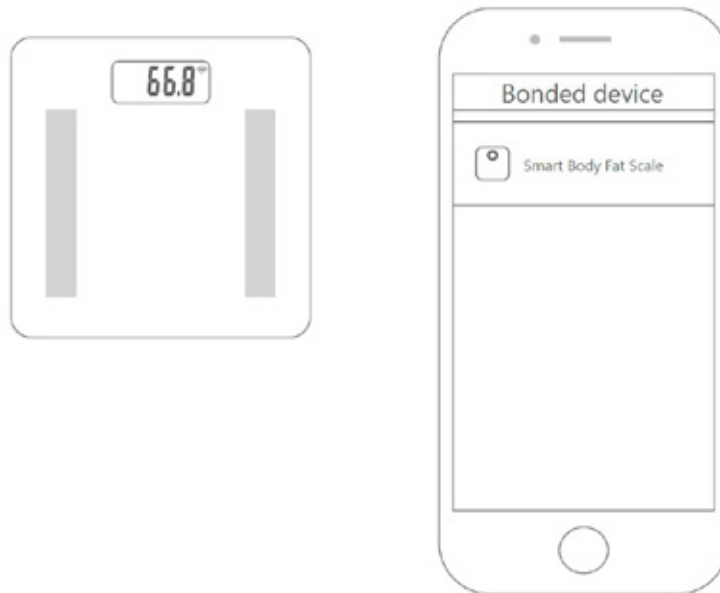
5. START MEASUREMENT

With the App, then put it on a solid and flat floor, and step on it for the measurement. Once the weight is stable and flashing, weigh measurement will be proceeded. You will review your weight and body fat display on the screen for several time, then the scale will be off automatically.


Do not put the scale on a carpet or other soft surface, it will lead to a wrong weight measurement.

- The first measurement has to be taken with Bluetooth connection, the first weight data will be reference for next auto user recognition.
- To make a proper measurement of body fat, please keep your bare feet on the electrodes of scale. The scale platform should be clean and dry.
- If you want to review all measurement data, open the App, and check them in detail interface.

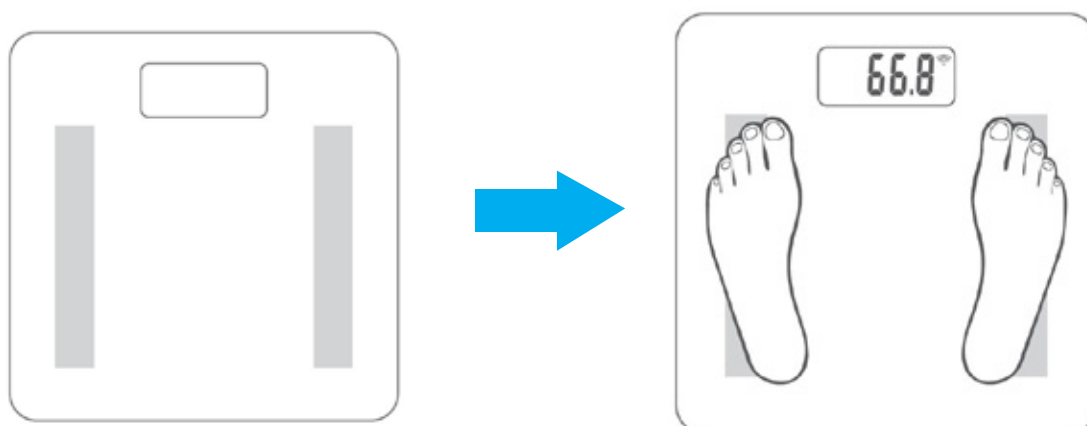6. MULTIPLE USERS SUPPORTED

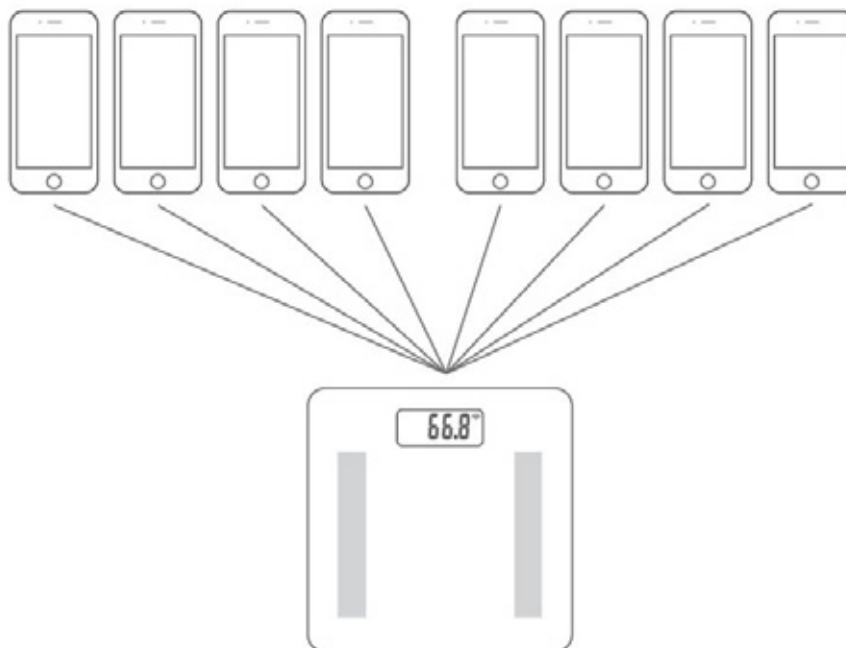8 accounts can be connected with a scale.

7. WARNING INDICATIONS

	<p>Fat Measuring</p> <p>Fat measurement is in process, please don't step off the scale, and make sure your bare feet is touching with the electrodes.</p>
	<p>Low Fat% Indication</p> <p>The Fat% is too low. Please lead a more nutritional diet and take good care of yourself.</p>
	<p>High Fat% Indication</p> <p>The Fat% is too high. Please watch your diet and do more exercise.</p>

	<p>Calibration</p> <p>Error exists, please re-test to get the correct result.</p>
	<p>Over-load Indication</p> <p>The weighing subject on the platform exceeds the maximum capability of scale. Please step off to avoid damage</p>
	<p>Low Battery Indication</p> <p>The battery power is running low, please replace with a new battery.</p>

8. ABOUT THE BODY FAT SCALE

The product is a smart body fat scale with a high-precision strain gauge sensor, which can help you to store and keep track of your weight on smart phone. This scale can transmit your weight data to your smart phone via Bluetooth technology, as well as make graph of your weight results, from which you are able to keep track of the trend anytime. Download the App on your phone to use Bluetooth scale, you will start your first step of weight control.

9. ADVISE FOR USE & CARE

1. Do not use the scale in humid, overheat environment. (Temperature range from 5°C to 35°C)
2. Stand on the scale steadily and touch the electrodes on the platform with your bare feet. Tip over will lead to a fall when you step on the edge of one side of the scale.
3. Do not drop, shock or strike the scale as it is a precise instrument.
4. Keep the scale clean. Do not use acid cleaning agent.
5. Do not use the scale on the place with intense shake. Keep standing firmly on scale when in use.
6. To ensure the accuracy of the data, it is recommended that each measurement be carried out at the same time of day.
7. Fat measurement is not available for children under 10 years of age and adults over 100 years of age.
8. When the battery runs out, please remove the battery in time to prevent leakage of the battery lead scale body corrosion;
9. The use of a fat measurement function is not permitted for pacemakers or other implanted medical devices.

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marca propiedad de | est une marque de | trade mark propriety of: